


CVX-RIPE[®]

Simultaneously provides gentle pressure at the internal and external os, as well as within the cervical canal.



- Conical shaped balloons, designed with no gap in between, form an hourglass shape when placed in the patient.
- Soft, silicone “nose” and catheter body combined with a pliable stylet to balance patient safety with ease of insertion.
- Catheter is packaged with a 20 mL syringe for clinician convenience.
- Luer-activated inflation ports for easy inflation.
- Not manufactured with Latex, DEHP or BPA
- Catheter tubing is 18 French.
- 5 year shelf life

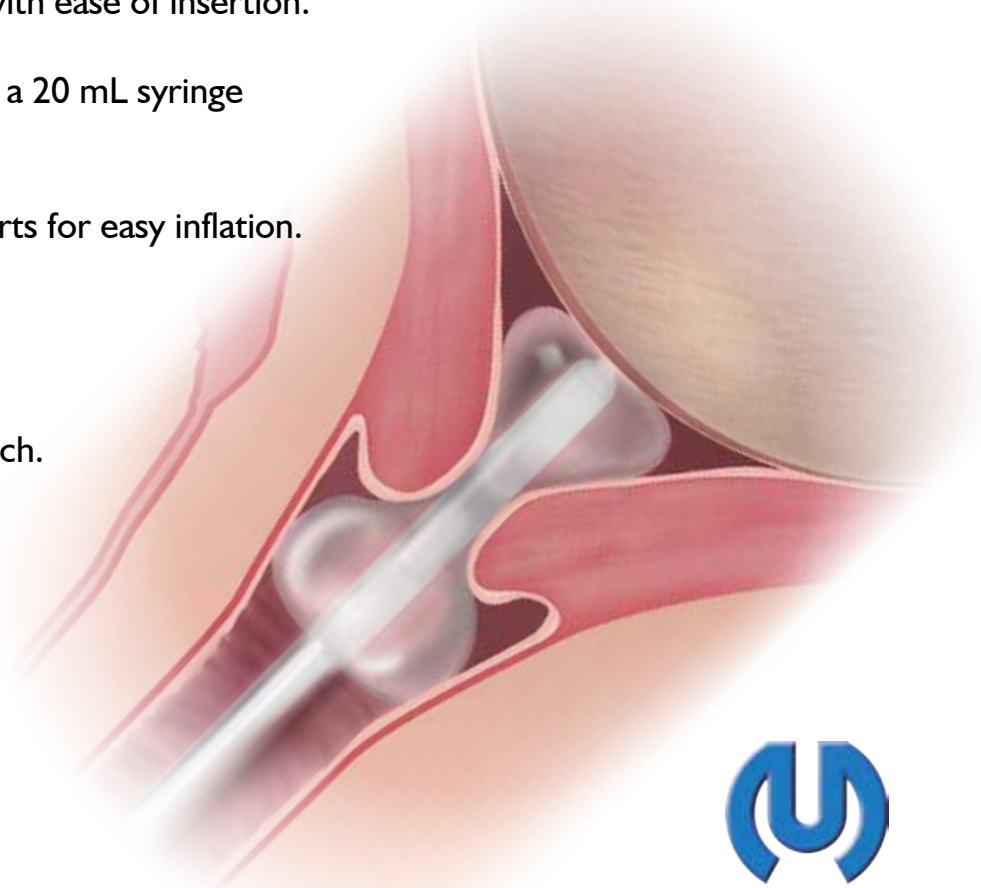

For more information or to place an order,
please contact Utah Medical Products, Inc at
Toll Free: **(800) 533-4984**
Web: **www.utahmed.com**

Part Number: CVX-100
Packaging: 10 per box



UTAH MEDICAL PRODUCTS
FEMCARE-NIKOMED

